

β -Catenin

Cat. # CM1181

Host Mouse Monoclonal IgG1


Size 100 μ l

Background:


β -Catenin is a 92 kDa protein that binds to the cytoplasmic tail of E-Cadherin. The cadherins, transmembrane adhesion molecules, are found with catenins at adherens junctions. Deletions in the cytoplasmic domain of E-Cadherin eliminate catenin binding and result in a loss of cell adhesion. Tyrosine phosphorylation of β -Catenin can regulate its interaction with critical components of adherens junctions. Both Fer and Fyn kinases phosphorylate tyrosine 142 *in vitro*. Overexpression of these kinases in epithelial cells disrupts interactions between α - and β -Catenins. The phosphorylation of tyrosine 142 may act as a switch from the transcriptional to the adhesive role of β -Catenin. Src family kinases can also phosphorylate tyrosine 86 and 654 in β -Catenin. The Tyr-654 phosphorylation regulates β -Catenin binding to E-cadherin. Thus, site-specific tyrosine phosphorylation of β -Catenin may regulate protein-protein interactions leading to changes in cell adhesion.

References

- Ozawa, M. et al. (1990) Proc. Natl. Acad. Sci. USA 87:4246.
 Roura, S. et al. (1999) J Biol Chem. 274(51) :36734.
 Piedra, J. et al. (2003) Mol. Cell. Biol. 23(7):2287-2297.
 Brembeck, F.H. et al. (2004) Genes Dev. 18(18):2225-2230.

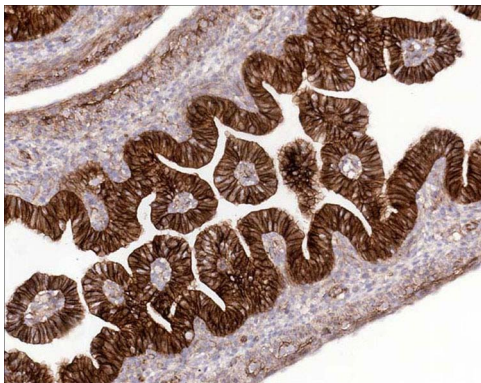

Formalin fixed, citric acid treated paraffin sections of embryonic Rat E16 intestines. Sections were probed with anti- β -Catenin (CM1181) then anti-mouse:HRP before detection using DAB. (Images provided by Carl Hobbs and Dr. Pat Doherty at Wolfson Centre for Age-Related Diseases, King's College London).

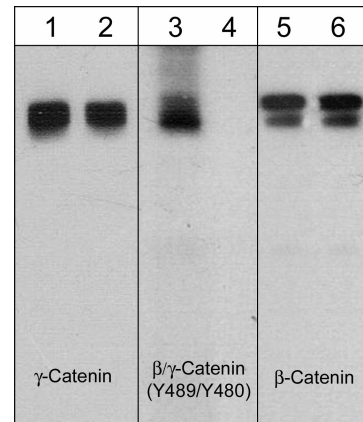

Western blot analysis of A431 cells stimulated with pervanadate (1 mM) for 30 min (lanes 1, 3, & 5) then treated with alkaline phosphatase (lanes 2, 4, & 6). The blot was probed with anti- γ -Catenin (CM1111), anti- β -Catenin (Tyr-489) conserved site (CP2961), or anti- β -Catenin (CM1181).

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

web: www.ecmbiosciences.com

telephone: 859-879-2075

email: info@ecmbiosciences.com

toll-free: 1-800-859-8202

β -Catenin

Cat. # CM1181**Host** Mouse Monoclonal IgG1**Size** 100 μ l**Immunogen:**

Clone (M118) was generated from a recombinant protein corresponding to amino acid residues in the C-terminal region of human β -Catenin. This peptide sequence is highly conserved in rat and mouse β -Catenin.

Buffer and Storage:

Mouse monoclonal antibody purified with protein A chromatography is supplied in 100 μ l phosphate-buffered saline, 50% glycerol, 1 mg/ml BSA, and 0.05% sodium azide. Store at -20° C. Do not aliquot. Stable for 1 year.

Applications:

WB	1:1000	ICC	1:250
ELISA	1:2000	IHC	1:500
IP	1:100		

End user should determine optimal dilution for their particular applications and experiments.
Western blot membranes were incubated with diluted antibody in 5% non-fat milk, PBS, 0.04% Tween20 for 1hour at room temperature.

Specificity:

The antibody detects a 92 kDa* protein corresponding to the molecular mass of β -Catenin on SDS-PAGE immunoblots of A431 and Hct116 src transformed cells. In addition, this antibody recognizes only β -Catenin in immunoprecipitations using anti- β -Catenin versus γ -Catenin.

*All molecular weights (MW) are confirmed by comparison to Bio-Rad Rainbow Markers and to western blot mobilities of known proteins with similar MW.

Related Products:

- CP1191 β -Catenin (Tyr-86), phospho-specific Rabbit Polyclonal
- CP1081 β -Catenin (Tyr-142), phospho-specific [Conserved site] Rabbit
- CP1201 β -Catenin / γ -Catenin (a.a. 649-661) Rabbit Polyclonal
- CP1061 β -Catenin (N-terminal) Rabbit Polyclonal
- CK6120 β -Catenin Phospho-Regulation Antibody Sampler Kit
- CK6230 δ 1-Catenin Phospho-Regulation Antibody Sampler Kit

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

web: www.ecmbiosciences.com

telephone: 859-879-2075

email: info@ecmbiosciences.com

toll-free: 1-800-859-8202

ECMBiosciences