

PDK1 (N-terminus)

Cat. # PP1411

Host Rabbit Polyclonal

Size 100 µl

Background:

3-Phosphoinositide-dependent kinase 1 (PDK1), also known as PKB kinase, was identified as the activator of the survival kinase Akt/PKB. Several important substrates of PDK1 include p70S6 kinase, PKAs, PKCs, SGKs, RSKs, and PAKs. PDK1 is a member of the AGC superfamily of serine/threonine kinases. Through the phosphorylation of downstream kinases, like Akt, PDK1 has been shown to be involved in several different cell functions, such as protein synthesis, cell survival, glucose metabolism, and cell adhesion and migration. The regulation of PDK1 occurs through lipid second messengers and phosphorylation. Multiple serine sites are phosphorylated on PDK1. Serine 241 phosphorylation is required for PDK1 activity, while serine 396 has been implicated in PDK1 nuclear translocation. Tyrosine phosphorylation may also regulate PDK1 activity. Tyrosines 9 and 373/376 are phosphorylated by c-Src in vitro. Tyr-373/Tyr-376 are important for PDK1 activity, while Tyr-9 phosphorylation permits Tyr-373/Tyr-376 phosphorylation by c-Src. In addition, Tyr-9 may be important during angiotensin-II-induced focal adhesion formation.

References

- Park, J. et al. (2001) J Biol. Chem. 276(40):37459.
Taniyama, Y. et al. (2003) Mol Cell. Biol. 23(22):8019.
Scheid, M.P. et al. (2005) Mol. Cell. Biol. 25(6):2347.

Immunogen:

PDK1 synthetic peptide (coupled to carrier protein) corresponding to amino acids in the N-terminus of human PDK1. This sequence is highly conserved in rat and mouse PDK1.


Applications:

WB 1:1000
ELISA 1:2000

End user should determine optimal dilution for their particular applications and experiments. Western blot membranes were incubated with diluted antibody in 5% non-fat milk, PBS, 0.04% Tween20 for 1 hour at room temperature.

Related Products:

PP1431 PDK1 (Tyr-9), phospho-specific Rabbit Polyclonal
PM1461 PDK1 (C-terminal Region) Mouse Monoclonal
PX1415 PDK1 (N-terminus) Peptide
PX1435 phospho-PDK1 (Tyr-9) Peptide

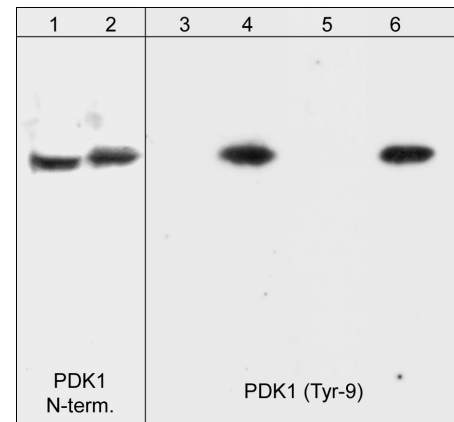

Western blot image of A431 cells untreated (lanes 1 and 3) or treated with pervanadate (lanes 2, 4, 5 & 6). Blots were probed with anti-PDK1 (PP1411) or anti-PDK1 (Tyr-9) (PP1431). The latter was used in the presence of no peptide (lane 4), phospho-PDK1 (Tyr-9) peptide (lane 5), or an unrelated phosphotyrosine peptide (lane 6).

Buffer and Storage:

Rabbit polyclonal, affinity-purified antibody is supplied in 100µl phosphate-buffered saline, 50% glycerol, 1 mg/ml BSA, and 0.05% sodium azide. Store at -20°C. Do not aliquot. Stable for 1 year.

Specificity:

This antibody was affinity purified using PDK1 peptide. The purified antibody detects a 65 kDa* protein corresponding to the apparent molecular mass of PDK1 on SDS-PAGE immunoblots of human A431 cells and mouse brain.

*All molecular weights (MW) are confirmed by comparison to Bio-Rad Rainbow Markers and to western blot mobilities of known proteins with similar MW.

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

www.ecmbiosciences.com
telephone: 859-879-2075
toll-free: 1-800-859-8202
tech: info@ecmbiosciences.com

ECMBiosciences

Rev 5/28/2008